



The Rebel is designed for users that find it difficult to maneuver a regular rollator or require extra support. The forearm supports are adjustable in height and angle find the optimal positioning for comfort and stability.

Outstanding Stability

The Rebel RA is the ideal rollator for pain sensitive users suffering from arthritis, Parkinson's or many other neurological disorders. It features separate forearm supports and large ball-bearing mounted wheels that ensure walking is easy, even on uneven or rough surfaces. The Rebel RA can be used both indoors and outdoors and features a wide assortment of accessories including a seat and basket¹.

Unique Handles

The Rebel RA features unique front-facing, angle adjustable handles that allow users to find the optimal handle position for comfort and stability. Adjustment is easy with

integrated measurements right in the handle tube. The ergonomic design of the grips ensures even distribution of weight. The Rebel RA can also be equipped with a one-hand brake to allow a user to control hand functions on either the left or right hand side.

Adjustable Forearm Supports

The Rebel RA features very unique forearm supports that are adjustable in both height and direction. This ensures that the forearm supports can be positioned in the best possible way for user comfort and stability.

¹Standard feature in Canada

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone : +61 1300 132 736
Fax : + 61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES


8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories/Features



Tray

Practical for resting lightweight items. Maximum load 2kg/4.4lb. Part no. 10601

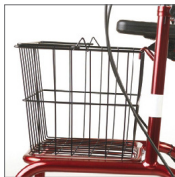


Basket

Included with rollator in Canada. Ideal for shopping. Part no. 10621



Slow-Down Brake

Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no.: 10620

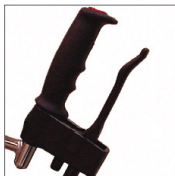

One-Hand Brake

Allows user to control hand functions with one hand (installed left or right). Special order item. Part no. 10653




Cane holder

Simplifies transport of a cane or umbrella. Part no. 10603



Reflector Kit

Includes 30 self-adhesive reflectors. Part no. 91090

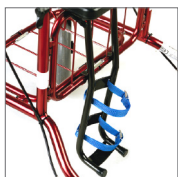

Oxygen Holder

Holds oxygen tube firmly in place. Part no. 10615



Curb Climber

Helps raise rollator onto a curb. Part no. 10612

This product is manufactured in accordance with all relevant ISO standards.



Technical Specifications

	Rebel RA Low 67	Rebel RA Low 72	Rebel RA Standard 67	Rebel RA Standard 72
Part no.	223111B	224111B	223211B	224211B
Handle Height (Min - Max)	86cm - 104cm / 33.8" - 40.9"	86cm - 104cm / 33.8" - 40.9"	104cm - 122cm / 40.9" - 48"	104cm - 122cm / 40.9" - 48"
Seat Height	54cm / 21.3"	54cm / 21.3"	62cm / 24.4"	62cm / 24.4"
Width between Handles	40.6cm / 16"	40.6cm / 16"	50.8cm / 20"	50.8cm / 20"
Wheel Size	20cm x 4cm / 7.9" x 1.6"	20cm x 4cm / 7.9" x 1.6"	20cm x 4cm / 7.9" x 1.6"	20cm x 4cm / 7.9" x 1.6"
Product Weight	11.3kg / 24.9lb	11.5kg / 25.3lb	11.4kg / 25.1lb	11.6kg / 25.6lb
Maximum User Weight	125kg/275.6lb	125kg/275.6lb	125kg/275.6lb	125kg/275.6lb

Follow Us on Social Media



facebook.com/humancarecda



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fjINh3zTjRxOP1IDf2Da0